


VisualEyes

User Manual

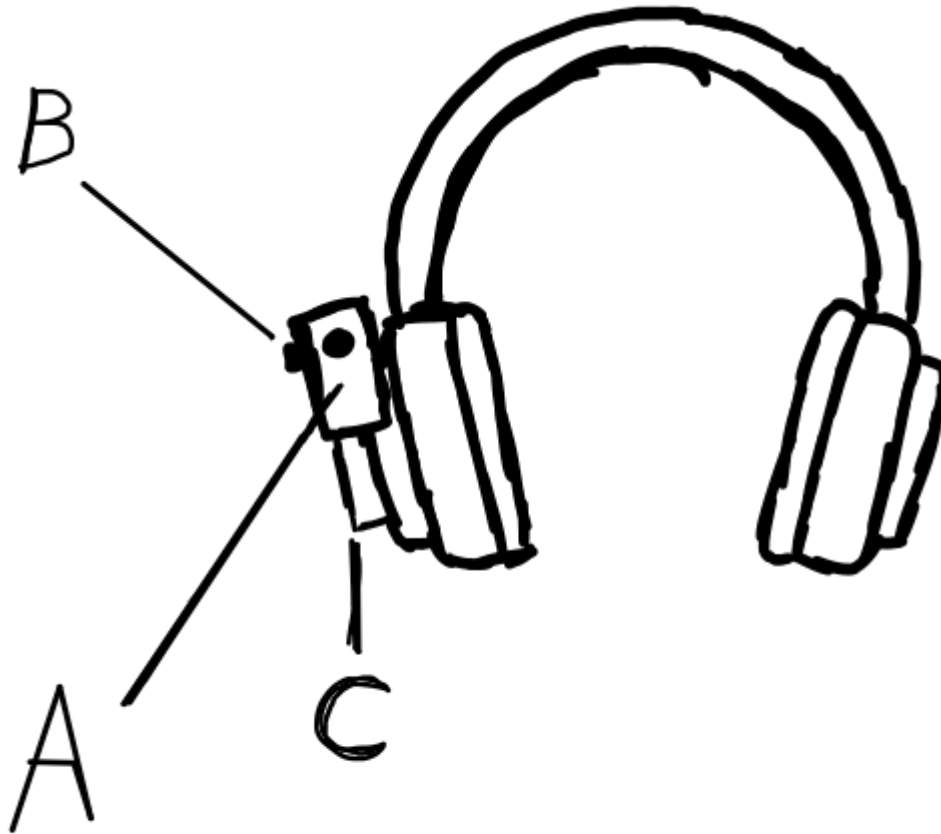

A - Camera

B - Action Button

C - Battery Casing

To use the VisualEyes you first place three AA batteries inside of the battery casing. Then you mount the casing onto a headwear of choice, such as a headset or hat. Then you press the camera until you feel a light click to turn the device on.



User Manual

After the VisualEyes is turned on, you can connect it to your phone through bluetooth. The device is named "Visual Eyes Camera XXX". When the device is connected, if you have the VisualEyes app on your phone opened, when you press the Action Button the phone speaks what the camera recognizes. Of course you can also connect your phone to an external speaker to have it speak through another device.

When the VisualEyes is not in use, remove the batteries from the casing to extend the life of the batteries and the camera.

If the device starts smoking or emitting a burning scent, remove the batteries from the casing or give a light tug on the wires to separate camera from the batteries. Failure to do so may lead to permanent damage to the camera or the product catching fire.



User Manual

Troubleshooting:

The camera is not listed under bluetooth devices:

- Ensure that the cables are connected
- Ensure that the batteries are present
- Press the camera until a light click is felt

VisualEyes

User Manual

VisualEyes App Manual

Once you open the app, the image below shows the default homescreen with certain buttons.

1. This button turns the Text-to-Speech module on. Once you click it, the icon will show a stripe through it, indicating it is turned off. When there is no stripe, it is on.
2. This box shows the output of what the AI recognizes.
3. The "Scene" button, marked with a camera icon next to it, is used to analyze a scene. Once you are pointing your camera at the object you want to analyze, click the "Scene" button. The AI module will detect the object, and the result will be displayed in the box above it.
4. The "Read Text" button allows you to detect text on different surfaces. Once you are pointing your camera at the surface you want to scan, click the button and the results will be displayed in the black box above it.

VisualEyes

User Manual

